

# Coatron® A6

Fully automated coagulation analyzer

- Measurement**
  - ▶ Clotting, chromogenic, immunological
  - ▶ Quarter test volume (75 µL)
  - ▶ High speed testing
  - ▶ Biphasic waveform analysis
  - ▶ Profile testing in batch or random access
  - ▶ Preset for 9 different profiles
  - ▶ Single or double determination
  - ▶ Automatic reflex testing
- Sample**
  - ▶ 3x 15 primary tubes
  - ▶ Positive sample identification (ID and rack position) with barcode
  - ▶ Cap Piercing for BD Vacutainer® or Sarstedt Monovette®
  - ▶ Micro tubes like Eppendorf® 2 mL for pediatric samples
  - ▶ Emergency STAT orders
  - ▶ Sample continuous loading in random access mode
- System**
  - ▶ 6 optic channels at 405 nm with auto light sensing and extraneous light elimination to reduce interferences such like bilirubin or sunlight
  - ▶ Integrated laser barcode reader for samples and reagent
  - ▶ Bi-directional Interface for PC-Link to LIS
  - ▶ Integrated graphical thermal printer
  - ▶ Heated and cooled positions for reagents, controls, buffers
  - ▶ Integrated drawers for consumables
  - ▶ System configuration as “OPEN” or “CLOSED”
  - ▶ Longlife Service Interval (every 100.000 tests)
- Quality Control**
  - ▶ Onboard colled positions for up to 3 control plasma
  - ▶ Auto QC flagging for erroneous results
  - ▶ One click report with Levey-Jennings graph and Westgard analysis
  - ▶ LOT. number (control + reagent) specific QC analysis is possible
- Software**
  - ▶ Intuitive TECO familiar graphical user interface
  - ▶ Multi language dialogue (En, Ger, Sp, It, China, Greek, Russia,...)
  - ▶ Free programmable test protocols to change or adapt new tests
  - ▶ Calculation of activity %, INR, Ratio, g/L, mg/dL and more
  - ▶ Calibration curve can be identified with up to 6 points per test.
  - ▶ Touchscreen Control (supported)
- Special Functions**
  - ▶ Automatic plasma dilutions up to 1:100
  - ▶ Automatic test calibration routines
  - ▶ Automatic test skip, if reagent run out
  - ▶ Database for 1 Mio patient history
  - ▶ Test calibration management for each reagent lot
  - ▶ Each result can be traced back to reagent lot and calibration data.
  - ▶ Patient monitoring for long time periods
  - ▶ Reaction curve monitoring
- Economic**
  - ▶ Quarter test volume to 75 % cost savings
  - ▶ Reagent dead volume is below 200 µL
  - ▶ Very economic consumption of rinse and cleaning solution
  - ▶ All critical system parameters are monitored and make the Instrument nearly free of service.

# Coatron® A6

Fully automated coagulation analyzer



- Global:**

PT  
aPTT  
Thrombintime  
Fibrinogen (Claus & PT derived)  
Heparin
- Technical specification**

Dimension  
50 x 95 x 80 cm

Weight (w/o packaging)  
52 kg

Power input  
90–240 Vac/50–60 Hz

Ambient conditions  
18 °C – 24 °C  
max. 60 % rel. humidity  
no direct sunlight
- Factors:**

Factor II  
Factor V  
Factor VII  
Factor VIII  
Factor IX  
Factor X  
Factor XI  
Factor XII

## Thrombophilia:

Protein C  
Protein S & free Protein-S  
Factor V Leiden (APC-R)  
Antithrombin  
D-Dimer & FDP  
LA Screen (DRVVT)  
LA Confirm (DRVVT)  
Plasminogen

## Calculated units

Activity %, Ratio, Polish ratio, INR (ISI system or INR calibration curve), g/L, mg/dL, ng/L, ng/mL, µg/L, IU/mL

## Optional available

- ▶ Enlarged waste management
- ▶ TECAM PRO LIS – Upgrade required for bidirectional ASTM interface to LIS

## Delivery package

- 1 Coatron A6
- 1 TECAM PRO CD + license
- 1 Power cord
- 100 Cuvettes, 6 pos/ea
- 10 Reagent container Ø 22,5
- 10 Reagent tubes Ø 13
- 10 Plasma tubes Ø 13
- 4 Magnetic stirrer
- 1 Rinse solution
- 1 Clean solution
- 1 Thermal paper roll
- 1 Waste container
- 2 Reagent adapter Ø 30.0
- 2 Reagent adapter Ø 24.5
- 2 Reagent adapter Ø 22.5
- 1 Operation manual
- 1 Download cable


TECO

Innovation in Coagulation

TECO Medical Instruments,  
Production + Trading GmbH

Dieselstraße 1  
D-84088 Neufahrn NB.  
Germany

Fon: +49-(0)87 73 - 707 80 0  
Fax: +49-(0)87 73 - 707 80 29  
Email: info@teco-gmbh.com  
Web: www.teco-gmbh.com

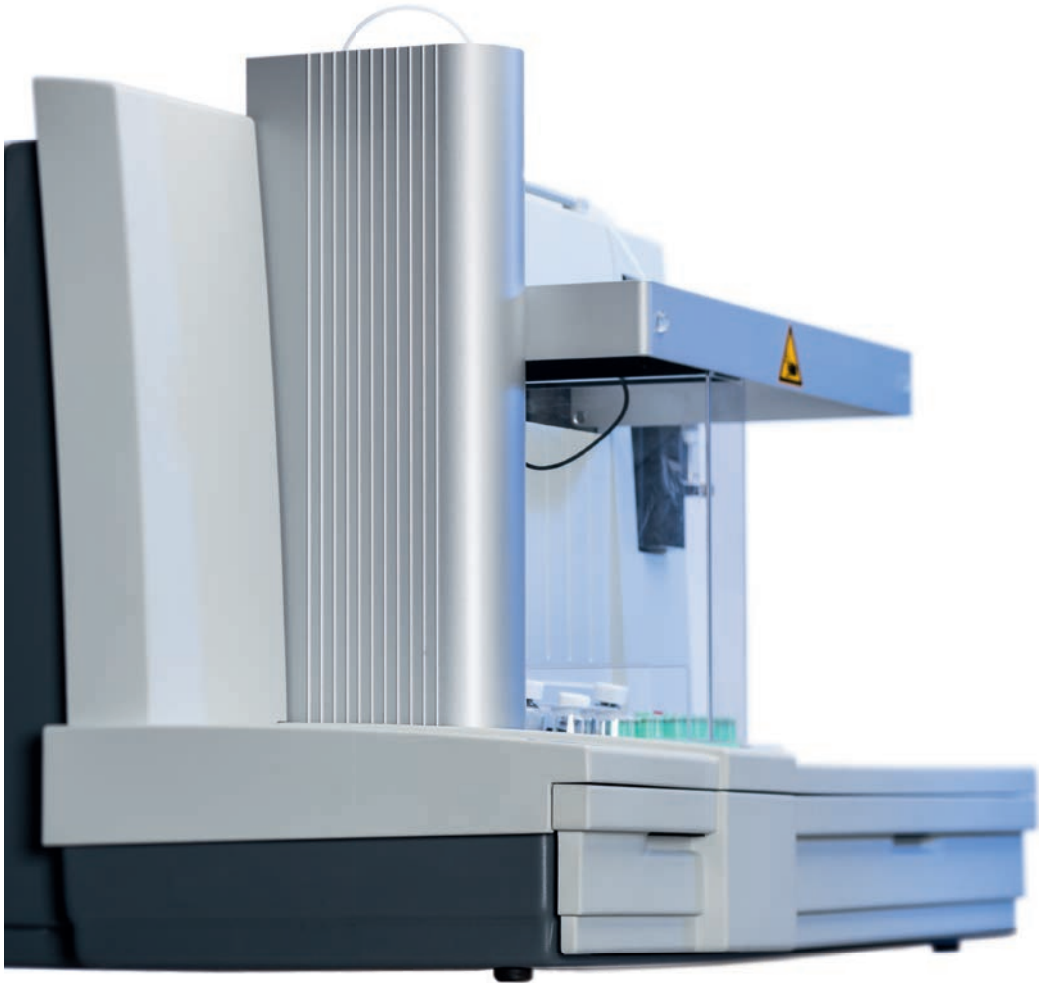

Made in Germany  
REF: 80 920 00

# Coatron® A6

Fully automated coagulation analyzer


- ▶ Clotting, chromogenic, immunological
- ▶ Up to 75 % lower running costs (quarter volume, low service)
- ▶ 45 patient, 34 reagent and 40 Cuvettes (240 tests)
- ▶ Integrated printer
- ▶ Microtubing for reagent and sample
- ▶ Integrated barcode reader
- ▶ Cap piercing
- ▶ STAT orders
- ▶ Patient continuous loading
- ▶ Reflex testing
- ▶ Automatic test calibration
- ▶ Automatic sample dilution
- ▶ Database for up to 1 Mio results (including reaction curve)
- ▶ Test calibration management
- ▶ Each result can be traced back to reagent lot and calibration
- ▶ QC reports with Levey-Jennings graph and Westgard analysis
- ▶ Multi languages
- ▶ Reagent barcode reading
- ▶ Reaction curve monitoring
- ▶ Programmable protocols to adapt new reagent or test
- ▶ Fully bidirectional LIS connectivity (RS232 or network)
- ▶ System monitoring for reagent, consumables and service
- ▶ Small footprint
- ▶ Longlife Service Interval (100.000 tests)
- ▶ System can be configured as “closed” or “opened”

TECO

Innovation in Coagulation

# Coatron® A6

Fully automated coagulation analyzer

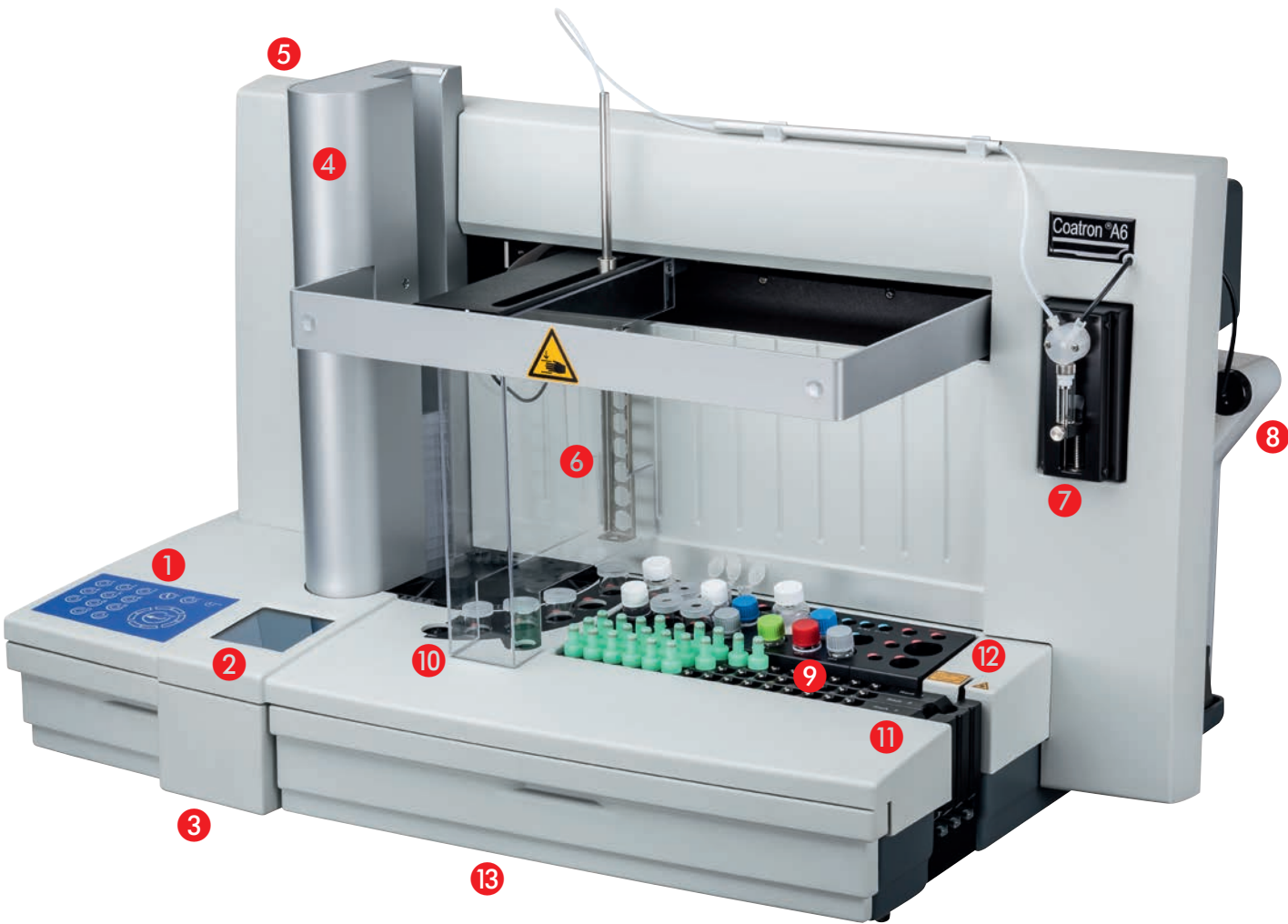

- 1 Keyboard
- 2 Grafic display  
Blue backlight
- 3 Cuvette waste drawer
- 4 Cuvette tower  
Enough for 240 tests (40 cuvettes)
- 5 Printer (back side)  
For printout of results, system analysis, test setups etc.
- 6 Liquid handling  
x-y-z high speed robotic system with Cap Piercing probe and level detection.


- 7 Pump/syringe  
Maximum volume: 2500 µl  
Minimum volume: 3 µl  
Long life syringe  
Easy to replace and maintain
- 8 Rinse container  
Contains 1,25 Liter Rinse solution, enough for appr. 1500 PT tests
- 9 Reagent block  
23 positions heated  
12 positions cooled
- 10 System block  
Wash & Waste position

- 11 Sample racks  
3 racks with 15 positions each, suitable for most common primary tubes. Every position have individual bar code for positive sample identification
- 12 Laser barcode scanner  
For patient, reagent and service barcodes
- 13 Liquid waste drawer  
For collection of liquid waste (optional: external tank system)

# TECAM

Software for TECO instruments

- Touchscreen Control
- Database for 1 Mio results
- Bidirectional LIS interface
- Easy worklist management
- Sample continuous loading
- STAT orders
- Automatic patient result and reports
- Visualize optical reaction curve
- Test calibration management
- Levey-Jennings and Westgard rules

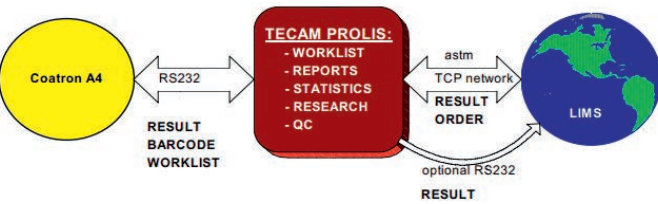

An ASTM interface is included, which allows to link with the laboratory information system (LIMS) over LAN network or RS232.

POS.: P02 ☒ STAT

Patient-ID: 2014000014169 ☐ Auto

Name (family, given): Mustermann, Paul ☐ Autodate

Date of birth: 10.01.1917 Age: 97 ☐ QC

Sex: m Set

Department: Oncology Set

Bed: 20

Test:

☒ PT ☐ DFIB ☒ APTT ☒ FIB

☐ AT ☐ TT ☐ D-Dimer ☐ HEP

☐ PC ☐ PS ☐ F5 ☐ F8

Define orders with a view clicks or by retrieving from LIS.

Coagulation Data Report  
Thrombosis & Hemostasis Research Center


|        |            |      |               |
|--------|------------|------|---------------|
| Name:  | Patient44  | SID: | 1069-4        |
| PID:   | 44         | Age: | 57            |
| Birth: | 19.06.1957 | Bed: | Oncology 2-15 |
| Sex:   | m          |      |               |

| Test            | Date             | Flag | Result    | Range       | Gratic. |
|-----------------|------------------|------|-----------|-------------|---------|
| PT (s)          | 23.09.2014 17:03 | A    | 21.7 s    | [10 - 16]   | ----->  |
| PT (%)          | 23.09.2014 17:03 | A    | 61.5 %    | [70 - 120]  | <-----  |
| PT (INR)        | 23.09.2014 17:03 | A    | 1.37 INR  | [0.8 - 1.2] | ----->  |
| APTT (s)        | 23.09.2014 17:03 | A    | 45.2 s    | [26 - 38]   | ----->  |
| D-Dimer (ng/mL) | 23.09.2014 17:06 | A    | 989 ng/mL | [0 - 150]   | ----->  |
| TT (s)          | 23.09.2014 17:06 |      | 19 s      | [16 - 24]   | ----->  |


29.10.2014 10:24

Receive results and create a patient report automatically.


TECAM software combines laboratory data management, quality control and research purpose in one. It connects the Coatron to the LIS. The worklist management allows defining orders much faster and flexible then on the instrument itself. Powerful filter functions will help to create any report. Levey-Jennings graph and Westgard analysis are used for quality control monitoring. Calibration data can be received from analyzer. Each result can be traced back to reagent lot and calibration.

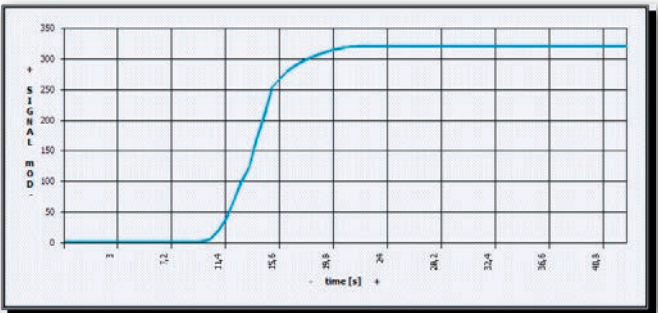

View reaction curve.

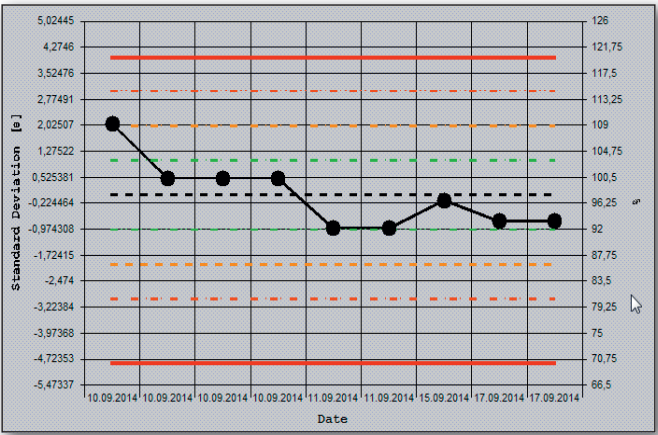

Create statistical reports with one click.

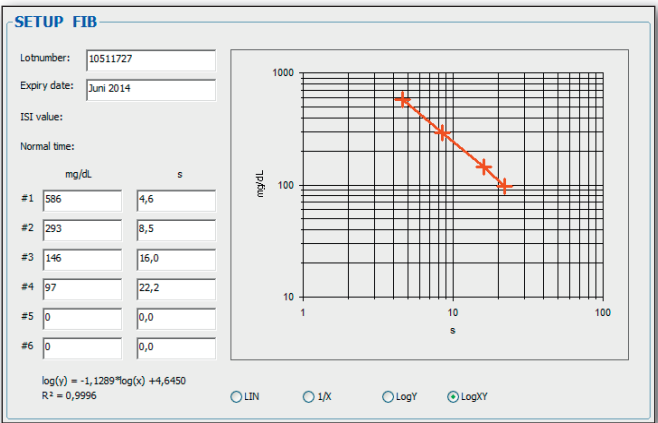

Each result can be traced back to calibration data.

# TECAM

Software for TECO instruments

- Recommended system requirements:  
CPU 1.5 GHZ, 2 GB RAM, Screen 1280 x 1024.  
Tested for Windows™ 7 or higher
- Instrument Interface  
Serial RS232 115200B, 8, 1, N  
Bidirectional for results including reaction curve and orders
- LIS Interface, unidirectional  
Serial RS232 115200B, 8, 1, N  
Forward results to LIS
- LIS Interface, ASTM  
TCP/IP network, LAN  
Bidirectional for forward results to LIS and request orders
- Intuitive handling  
One click principle for data collection, filtering, reporting or QC
- Touchscreen Control  
Operate with keyboard, mouse or/and touchscreen
- Multi lingual  
English, Spanish, Italian, Portuguese, Russian, Greek, Chinese, further languages possible
- Database  
Manage over 1 Mio results in one database. Powerful SQL filters will help to find and select any results.
- Traceability  
Each result can be traced back immediately to reagent lot and calibration data.
- Calibration data  
Test calibration can be import, export, edit for each method.
- Reporting  
Make cumulative reports by one click. Use SQL filters to generate specific reports (e.g. results of patient of all PT within the last 30 days).

- Order management  
Create, load, save, repeat, receive and send orders. Positive patient identification by barcode validation. USB barcode scanner supported. Up to 9 different test profiles. Preset for up to 6 different QC samples.
- STAT orders  
Each order can be marked as STAT and will be executed as soon as possible.
- Flexible multi filter  
Combine every data field like  
PID=1234 + Date=today + Test=FIB.
- Filter Presets  
Today; recent x-days; QC data; etc.
- Statistical analyse for QC or patients  
Filter any data (e.g. Patient XYZ, only PT, last 30 days), and visualize them with Levey Jenning diagram with one click. Analyze data with Westgard Multirules. Define QC or Patient ranges.
- Reaction curve analysis  
Each result is received together with optical reaction curve.
- Research  
Export all data, grids, charts to PDF, XLS, TXT or HTML.
- Security  
Operator can login as administrator or user in order to control and manage the database with different access rights
- Online Update